



www.gmedi.com



I.V FLOW REGULATOR




www.gmedi.com



www.gmedi.com

Product Features

- Modernism hands & one hand type free design
- Pressure sensitivity accuracy regulation membrane fully adjusting 0 ~ 250ml / h
- Stability of flow rate over the 24hours of infusion allows
- Soft tubing
- PVC FREE / DEHP FREE
- I.V Flow Regulator Protector (Option)
- I.V Spike & Dextrose Fixer (Option)



Option

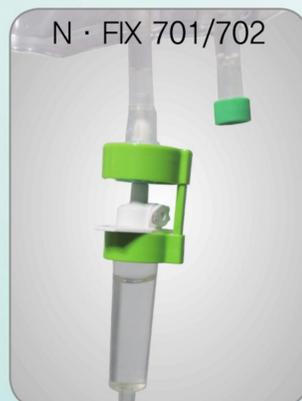

특허 제 10-1519304호

Filter



Round Type



Rectangular Type

(Single/Double
5 μ m, 1.2 μ m, 0.2 μ m)



Manufacture Address:

#401, Chunkang Apartment Shape Factory 218, Gurojungang-ro, Guro-gu, Seoul, 152-887, KOREA
TEL: +82-2-2632-9256 FAX: +82-2-2631-2659 E-mail: gmedir@yahoo.co.kr
http://www.gmedi.com http://eblog.gobizkorea.com/gmedi http://kr109859445.trustpass.alibaba.com/

* The above specifications are subject to change without prior notice.



www.gmedi.com

Medical Disposable 3 Way Stopcocks & Extensions


DEHP FREE / PVC FREE

DIEL™


www.gmedi.com



• Filter

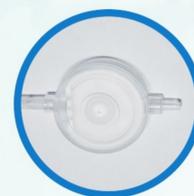

Round Type

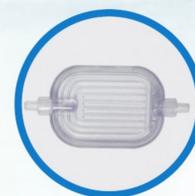

Rectangular Type
(Single/Double)
(5µm, 1.2µm, 0.2µm)

3Way stop cocks



1. CE and ISO certified
2. Made from polycarbonate for additional channels of administering I.V. fluids
3. Tap turns every 360 degree without limitation
4. Multiple channels for multiple infusion therapy
5. Designed to withstand pressure up to 5 bars

Part materials:

Part Name	Material
Threaded stopper	High density polyethylene
Stopcock tap	High density polyethylene
Tube	PVC free
Male rotator	ABS
Mediflexe body	Poly carbonate
Plain stopper	Poly propylene

3Way Extensions



Manufacture Address: